

Zemits Ultranexx Lite




| Introducing non-surgical
high intensity focused ultrasound
technology for smas-lifting

Technology


Zemits Ultranexx Pro is the latest double-feature platform for skin lifting and targeted fat reduction both for the face and the body. The device is known for its innovative technology that eliminates the need for surgery, does not cause injuries and does not leave scars.

Thanks to the Zemits Ultranexx Pro's faster-shot speed, the treatment is comfortable and results in no side effects, downtime, and scarring for most people. HIFU can be used to treat any area of the face and body with seven different cartridge depth options 1.5mm – 16mm.

Deliver long-lasting results to your clients and receive in-house training and lifetime support from Zemits.

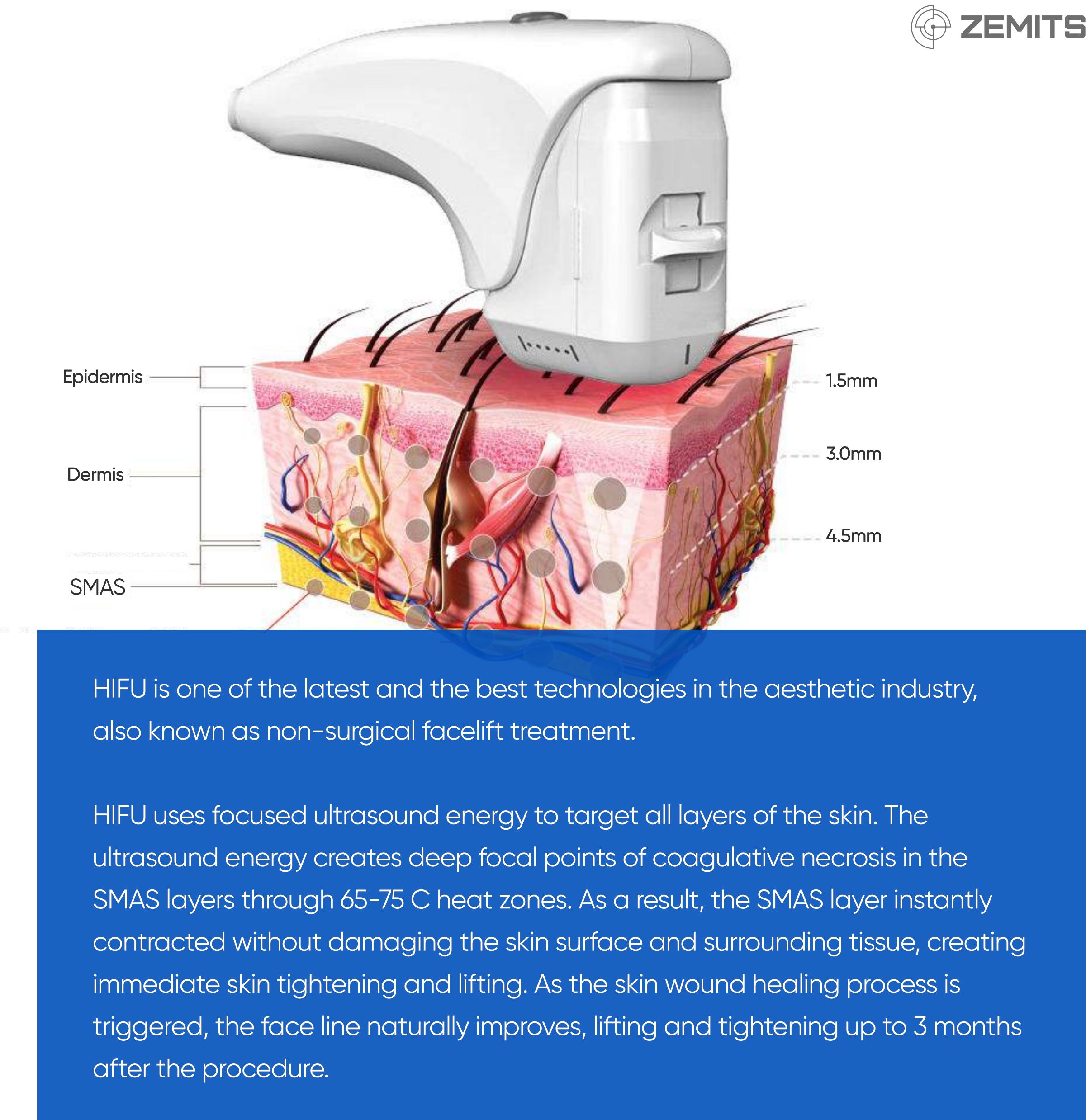


KEY ADVANTAGES

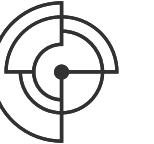-  Double-application system for both face and body treatments
-  The safest and most effective non-surgical technology for skin tightening and lift
-  Only one treatment is required for a lasting lifting effect
-  Easy-to-operate with SMART touch interface
-  Low consumables expenses
-  LOW Investment and HIGHLY profitable treatment
-  SMAS lifting treatment increased demand to expand clientele



Zemits Ultranexx Lite is recommended for the treatment of

- Increased discoloration of the skin;
- Decreased skin density;
- drooping skin in the cheek area
- Marionette lines
- Deep nasolabial folds
- Ptosis
- Double chin (flabby and sagging skin under chin);
- Mummy Tummy
- Localised fat pocket(s) which will not budge even with an active and healthy lifestyle
- Post weight loss (over stretched skin on the arms, legs, above the knees, abdomen, face, neck, décolleté)
- In combination with other body sculpting treatments

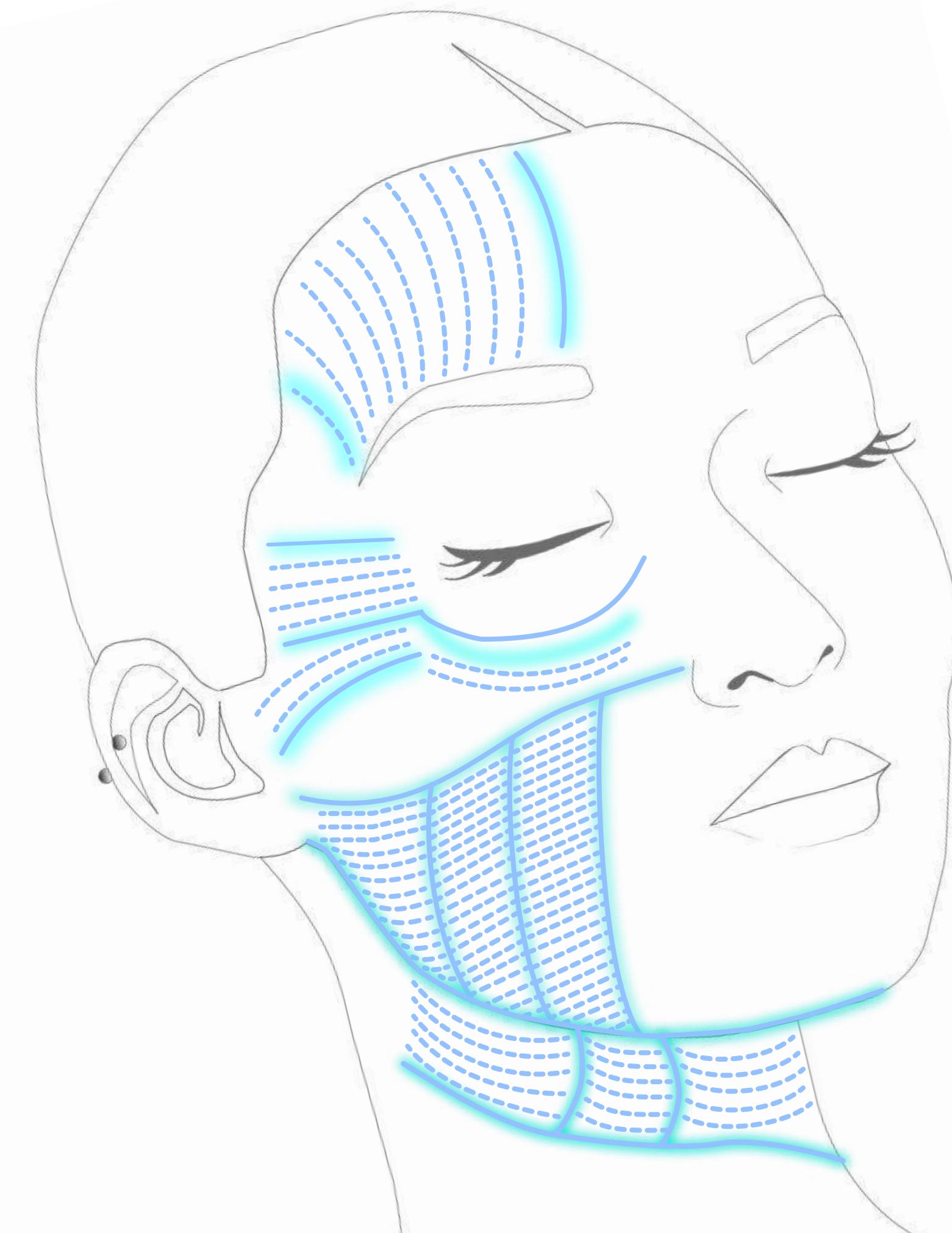

Zemits Ultranexx Lite Double-application system:

FACE AND NECK LIFTING


- Skin rejuvenation
- Reduction of wrinkles
- Tightening of sagging skin around the neck
- Lifting cheeks, eyebrows, and eyelids
- Improving the contours of the jawline
- Tightening of the décolletage
- One-off treatment
- Final result after 2-3 months and lasting up to two years



Zemits Ultranexx Lite Double-application system:

BODY SCULPTING

- Targeted fat reduction
- Reduce the appearance of cellulite
- Firm and tighten the skin
- Increase elasticity in the areas targeted
- Treatment zones: abdomen, thighs, flanks, hips and buttocks
- Safe alternative for liposuction surgery
- Visible results after 2 months post treatment

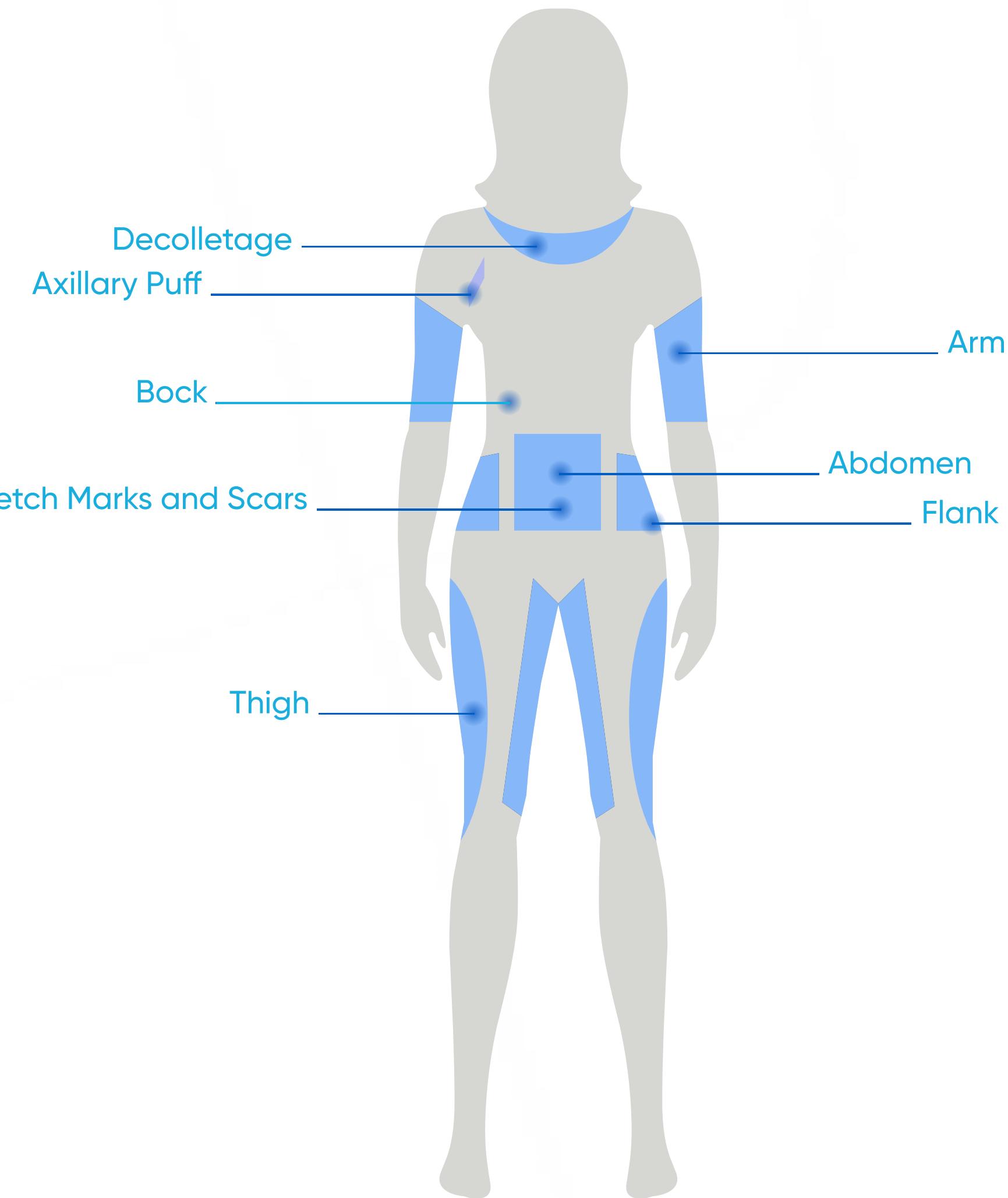

Treatments

EYES PROTOCOL (30 min)

This treatment helps to reduces puffiness, dark circles, and firm the muscles around the eye area. Start with Vacuum lymphatic drainage massage using Zemits Dazzle Skin oil. Remove all oil residue and prepare gel. Select the lowest depth cartridge depth for treatment on the eye area. Once both eyes have been treated, remove all gel from skin and apply a soothing eye mask. We recommend also using LED to promote the healing and rejuvenation.

MUMMY TUMMY PROTOCOL (45 min)

This treatment is designed to tighten, restructure and sculpt the body post pregnancy. To further enhance results we recommend starting with a manual/mechanical lymphatic drainage massage. Remove any oil residue, cleanse and mechanically exfoliate area lightly. Select the corresponding protocol cartridge depth for the treatment area, once the first pass is complete, follow with the corresponding cartridge to "layer" the following targeted depth. Once all areas have been treated, remove all gel from skin and apply an Intense Hydration moisturiser. We recommend inviting the client back for a course of pressotherapy treatments to further enhance the results.



Treatments

ANTI-AGEING PROTOCOL (90 min, Full face, neck & Decollete)

This treatment is designed to address signs of ageing using a combination of targeted depths. You can deep cleanse with the Ultrasonic Scrubber skin spatula first, removing any debris from the skin. Select the corresponding protocol cartridge depth for treatment area, once the first pass is complete, follow with the corresponding cartridge to "layer" the following targeted depth. Once all areas have been treated, remove all gel from skin and apply an Intense Hydration mask. We recommend also using LED to promote the healing and rejuvenation.

SCAR & STRECH MARKS PROTOCOL (30 min, targeted treatment)

Designed to soften and refine the appearance of scarring through targeted resurfacing. Use a Skin Scrubber Ultrasonic skin spatula for initial surface peeling. Select the corresponding protocol cartridge depth for treatment area, once the first pass is complete, follow with the corresponding cartridge to "layer" the following targeted depth. Once all areas have been treated, remove all gel from skin and apply an Intense Hydration mask. Once completed we recommend also using LED to promote healing and rejuvenation



INCREASE YOUR REVENUE AND RETURN YOUR INVESTMENT IN NO TIME
with new treatments that appeal to every client.

By investing in Zemits Ultranexx Lite you are investing in your future successful aesthetic business and getting a solid partnership with a multinational family you can trust to deliver quality, trustworthiness and performance. Zemits Ultranexx Lite tremendous earning potential without compromising technology as an affordable addition to any clinic or salon.

Offering Zemits facial treatments within your business will enable you to attract new clients, increase retention and significantly boost your revenue.

The potential to earn £250 000 in the next 6 months

based on two treatments per day

	Time for treatment	Packages	DAY	1 WK	1 MTH	6 MTHS
Full Face, Neck & Decollete Treatment	90 min	1 session	£2085	£10 425	£41 700	£250 200
Small area (Body)Treatment	30 min	3-6 session	£1185	£5925	£23 700	£142 200
Eye treatment	30 min	1 session	£360	£1800	£7200	£43 200
Mummy Tummy treatment	45 min	3-6 session	£1785	£8925	£42 840	£257 040



Types of cartridges for operation

With Zemits Ultranexx Pro you can use cartridges of different penetration depths.

With your device purchase you receive the 3 standard cartridges for facial treatments:

- 1.5mm
- 3.0mm
- 4.5mm

The deeper body cartridges 6.0 mm, 8.0 mm, 10.0 mm, 13.0 mm and 16.0 mm are also available for purchase.



1.5 mm



3.0 mm



8.0 mm



10.0 mm



4.5 mm



6.0 mm



8.0 mm



10.0 mm



13.0 mm



16.0 mm

Why Zemits?

Zemits UK - an American brand born in sunny California. A team of like-minded people seeking to provide professionals worldwide with aesthetic equipment and products to help every beauty business succeed.

The core of the Zemits brand consists of outstanding user-friendly equipment and ongoing marketing and client support.

With every product we offer:

- ⌚ 2 years warranty on the device
- ⌚ Free Training
- ⌚ Free Treatment Protocols
- ⌚ Marketing support materials
- ⌚ Finance options



Professional training with
Zemits qualified experts



Marketing Materials for
boosting your business



2-Year
Warranty

We at Zemits believe – time is money!
Although you are always welcome in our showroom, we understand that it's not always possible to find time to visit us.

So, trying to be flexible for your comfort, we provide a face-to-face VIRTUAL CONSULTATION session so you can learn about all the advantages of using Zemits advanced devices for your business!

Our experienced aestheticians will provide you with a personalised training plan for every piece of Zemits equipment.

Book your Virtual Demo face-to-face to get fully equipped with the confidence to start with Zemits systems in your business, regardless of location.

[CHECK REVIEW](#)[BOOK CONSULTATION](#)

OUR CONTACT DETAILS

We are delighted to invite you to our showroom, where you can receive a full consultation and try out the device.

Our Business Hours: 09:00 - 18:00 (GMT) Mon-Fri,

+44 (20) 45325922

Place an Order Online 24/7

zemits.co.uk

